

Service Instructions

SP014/03/I

Title: Cover at head end

Reason for update:

☐ Reliability☒ Infrequent ServiceUrgency: ☐ Immediate ☒ Within 12 monthsUpdate materials required? ☐ Yes ☒ NoMaterials free of charge? ☐ Yes ☒ NoReturn of parts? ☐ Yes ☒ No

Estimated completion time: 0,5 Number of CSE's:1

Customer application training? ☐ Yes ☒ No**Systems/Products affected/System identifying IVK**

Name	Material No.	Serial No.
Uroskop ACCESS LA without cassette	57 56 122	1001 - 1029 1031 - 1037
Uroskop ACCESS RA without cassette	57 56 130	2001 - 2013

Components affected/to be modified

Name	Material No.	Serial No.	Component status Affected
n.a.	n.a.	n.a.	n.a.

Remark:

Chg. Rev. No.: 115440
Name: D.Fürstenberg
Dept.: CS SD 24

© Siemens AG 2003

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

Chapter	Page	Revision
all	all	01

Document Revision Level

This document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

Disclaimer

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such service.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

1	General Information	1 - 1
	Uroskop ACCESS	1 - 1
	Reason for the Update	1 - 1
	Prerequisites	1 - 1
	Special Tools / Documents	1 - 1
	Ordering Information	1 - 2
	Return of Parts	1 - 2
2	Work Steps	2 - 1
	Work Steps	2 - 1
	Final Check	2 - 2
	Customer Information	2 - 2
	Final Work Steps	2 - 2
3	Changes to previous Version	3 - 1
4	Update Completion Form	4 - 1
	Completion Protocol / Update Completion Form	4 - 1

This page intentionally left blank.

Uroskop ACCESS

Uroskop ACCESS systems without cassette if liquid drops inside the table at head end.

Reason for the Update

Uroskop without cassette may have not the correct cover holder at table head end installed and liquid may drop inside the table.

Prerequisites

Order a new cover holder.

Special Tools / Documents

n.a.

Ordering Information

The following spare part has to be ordered from CSML:

for

Uroskop ACCESS LA without cassette no. 57 56 122

Cover holder **73 77 141**

for

Uroskop ACCESS RA without cassette no. 57 56 130

Cover holder **73 77 828**

Return of Parts

n.a.

Work Steps

- Open the cover at the table head end, see Fig. 1 (example for Uroskop ACCESS left hand version).



Fig. 1

- Check, if liquid is inside the table at table head, see arrows Fig. 1



Fig. 2

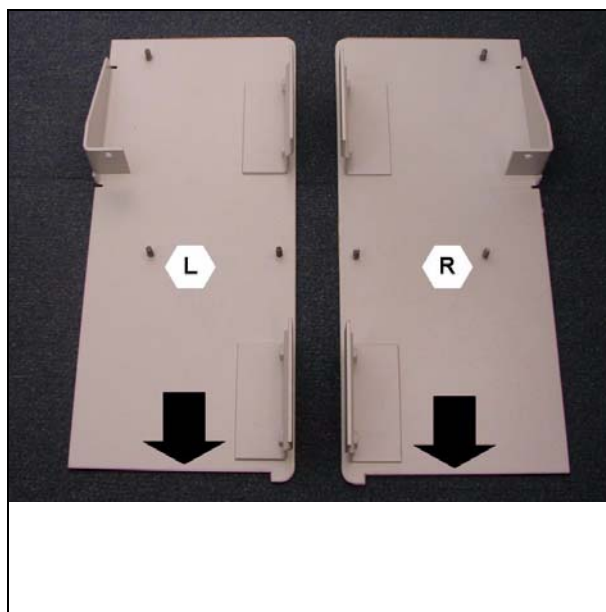


Fig. 3

- Check, if the cover holder is like Fig. 2
 - the edge at the holder have to be plain, see arrows Fig. 3.
 - "L" for Uroskop ACCESS LA, "R" for Uroskop ACCESS RA
- If needed, change the holder.
- Install the cover.

Final Check

- After assembly make sure that the rubber lip seals over the whole length of the cover.

Customer Information

n.a.

Final Work Steps

n.a.

No changes; initial publication.

This page intentionally left blank.

Completion Protocol / Update Completion Form

The update with the number **SP014/03/I** has been completed.

Material number:

Serial number:

Customer: Functional Location:

Customer No.:

Name (CSE): Telephone:

Country: Location:

Date: Signature:

Remark:

.....

NOTE

After completing the update, make a copy of this page, fill it out and file it in the corresponding System Binder/User Handbook.

This page intentionally left blank.